

# Golvo



# The complete mobile lift

Based on a world-patented lifting principle, Golvo is unique in terms of flexibility and versatility.

## Mechanical emergency lowering

Patented emergency-lowering function for quick and easy lowering.

## Handcontrol

Easy-to-use handcontrol with magnetic mount.


## Built-in charger

The built-in charger always stays in place on the lift. Simply charge the batteries by plugging the charger cable into a wall socket.

## Electrical emergency lowering and raising

In the event that the handcontrol for some reason stops working.



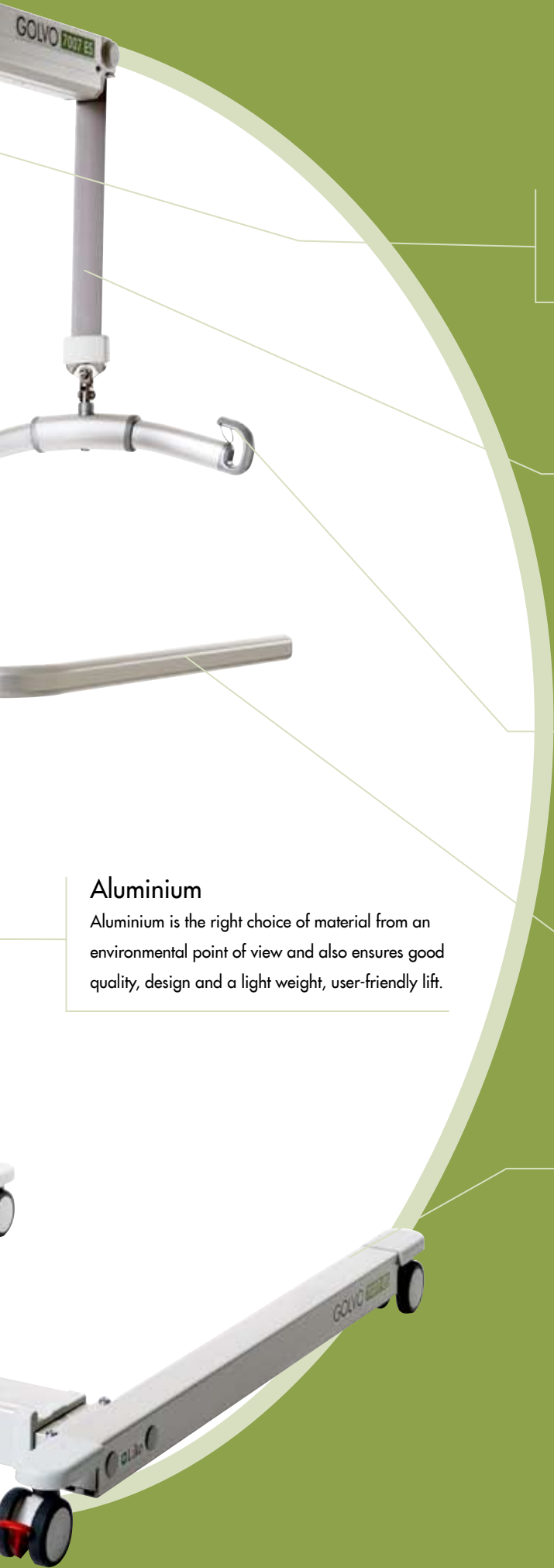

### Vertical lifting

The easiest and most space efficient lifting motion.

### Flexible lift strap

Adapt the lift to suit the person being lifted.

### Slingbar with safety latch

Prevents the slingbar from unhooking unintentionally, and makes the lift safe and secure.

### Retractable arm rests

For both the patient and the caregiver – added reassurance.

### Parallel base-width adjustment

Golvo's base can be widened and narrowed with the unique parallel adjustment, enabling access around beds or wide chairs etc.

### Aluminium

Aluminium is the right choice of material from an environmental point of view and also ensures good quality, design and a light weight, user-friendly lift.

# In your situation



## The "chameleon" of patient lifts

Unique in design and function, Golvo has been developed with an insight into the great variety of lifting needs that arise in the care and healthcare sectors. Golvo's flexibility and versatility, together with Liko's comprehensive range of functional slings and accessories, make it the ideal lift for most situations.

### ↑ Safety and rehabilitation

Together with Liko MasterVest or Liko Lift Pants, Golvo offers a popular means of risk-free gait training. Armrests give support and reassurance. MasterVest can also be used for visits to the toilet.

### ← Lifting children

Golvo's wide lifting range and versatility make it the ideal aid for lifting and stimulating children and young people. It can lift to high tables and lower right down to the floor. A wide range of lifting accessories specially adapted for children are available.

### ← Parallel base-width adjustment

Golvo's unique base can be widened and narrowed with the unique parallel adjustment. This allows the patient to be lifted into a wheelchair from the perpendicular and still achieve the correct sitting position.

### ← Lift with ease from the floor

Lifting a patient from the floor is surprisingly easy. The open base provides plenty of room, and the flexible slingbar suspension with one strap enables you to adapt the lift to suit each patient. The Liko Original High Back Sling gives the patient the best support.





Life is fantastic!

Everyone should be given the right to live life  
in the best way possible

# Easy getting to the toilet



## When the need arises

Lifting in the toilet area is often difficult and challenging. Golvo is versatile, even in hygiene areas where space is limited. There are slings for every situation. As always, the most appropriate sling is selected on the basis of the patient's condition and clothing.



### ↑ With Liko Original Sling (Mod. 10)

Often, toileting can be managed easily with one of Liko's basic sling models, e.g., Original or Universal.

The patient can be transferred to the toilet with a mobile hygiene chair.

### ← With Liko Hygiene Sling (Mod. 40, 41, 45)

Patients with good upper body strength can use Liko Hygiene Sling. Liko Hygiene Sling has a generous opening around the bottom that facilitates undressing during the lift.

### With Liko Hygiene Vest (Mod. 50, 55) →

Liko's unique Hygiene Vest solves most toileting problems. It provides unique advantages in terms of reassurance to both patient and caregiver. Many patients appreciate its stable, upright sitting position.

Hygiene Vests are available with head supports.

### ← Liko MasterVest (Mod. 60, 64)

Golvo's high lifting height makes it possible to raise a patient to the standing position. This facilitates undressing and dressing whilst toileting. Anyone capable of supporting their weight while standing can be lifted to a standing position with Golvo and Liko MasterVest.



# Respect and integrity

are important  
when you are lifting  
another  
human being





Liko OctoStretch

# Safe, secure lifts require more than just technology



Stretcher with Stretch Leveller



LikoScale with Quick-release adapter



# Liko Stretchers

## Rely on Golvo in every situation

Lifting a person from a supine or prone position places high demands on the lift and lifting accessories. Golvo's stability is reassuring. The lift is achieved safely and securely, mainly due to the unique parallel adjustment, plus the lifting motion is vertical, keeping the centre of gravity in the middle of the lift. Many lifting accessories are available for Golvo and the special needs of the acute ward. For example:

### ← Liko OctoStretch

OctoStretch, pictured here with a full lifting sheet, is an easy-to-use and much-appreciated aid in operating theatres and intensive-care units. Several types of lifting sheets are available for OctoStretch.

### LikoStretch 600 IC →

A horizontal lift for the most difficult and delicate patients. The lift strap is passed in under the patient with a pliable spatula, avoiding the need to turn the patient. Liko's unique strap-lock secures the patient in place, so there is no need to alter his or her position during lifting.

LikoStretch 600 IC with full sheet 1800 is appropriate for lifting and transferring patients who can be turned on their side.

### ← LikoScale

A user-friendly scale is available as an optional accessory for use with Liko's lifts. Weighing patients in the lift is practical and convenient. The patient sits comfortably in the sling whilst the caregiver reads the weight, clearly visible on a display at eye level. Mounted with the quick-release adapter, the scale can be easily moved when necessary.

### ← Liko Stretch Leveller

With the Stretch Leveller mounted on the horizontal lift, the balance can be adjusted quickly and easy, avoiding tipping of the patient.

### FlexoStretch

Designed for the heaviest lifting situations with the patient in a horizontal position. Lift sheet or lift straps are optional accessories. A Stretch Leveller is part of FlexoStretch.



LikoStretch 600 IC



FlexoStretch



# Down to the smallest detail

Golvo is the most versatile mobile lift on the market. Together with Liko's lifting accessories, Golvo is ideal in most lifting situations.



## The mobile overhead lift

Golvo is using the principles of an overhead lift. This gives unique flexibility.



## Vertical lifting motion

Gives maximum stability, since the centre of gravity is always directly in the middle of the lift.



## Parallel base-width adjustment

Golvo's base width has a unique parallel adjustment. This facilitates access in many different lifting environments, to and from the bed, toilet, x-ray table, etc.



## Adjustable lifting range

When lifting from the floor, the slingbar can be altered to the desired level. Golvo can reach all the way.



## Armrests

The retractable armrests give reassurance and stability for the user when gait training, toileting etc.



## Mechanical emergency lowering

Patented emergency-lowering function for quick, trouble-free lowering.

## Smart accessories

### Quick-release hook and Q-link

Liko has developed a system that enables quick and easy exchange of lifting accessories - without the need for tools. Change from a slingbar to a stretcher, or connect a scale in one simple operation. Pictured, LikoScale with quick-release hook TDM and Q-link on the lift strap.



### Extra charger station

With a wall-mounted charger station and an extra battery, reserve power is always available.

# Intelligent electronics

## Handcontrol

Two speeds for even control and base width adjustment.

## Built-in charger

The batteries and charger are contained in the one control box attached to the lift. Just plug the charger cable into a wall socket. Overcharge limiter terminates charging when the batteries are fully charged.

## Display for charge status

A clear audible and visual indication on the display tells you when the batteries need to be charged.



## Electrical emergency lowering and raising

If the handcontrol for some reason stops working, Golvo can be operated directly from the control box. The raising function allows you to make urgent lifts.

## Electrical emergency stop

Cut the power supply quickly by simply pressing a button.



## Technical data

Lifting capacity:	200 kg
Lifting speed:	Two speeds; 4.8 cm/s and 3.2 cm/s, respectively, without load
Batteries:	2 x 12V 2.9 Ah gel-type batteries (valve regulated rechargeable leadacid batteries)
Battery charger:	Built-in charger for 100-240 V, 50-60 Hz, max 400 mA
Material:	Anodized aluminium



## Safety latch

A simple, ingenious device that prevents the sling from unhooking unintentionally.



## Colours and sizes

Colour coded sling size chart added to co-ordinate with the Liko range of slings.



Golvo 7000 ES/7007 ES is a Class 1 medical-technical product that complies with the medical-technical directive MDD 93/42/EEC and the requirements of EN ISO 60601-1.

Made in Sweden

## Lifting accessories

Together with Liko's comprehensive range of accessories, Golvo offers solutions to virtually all lifting problems in the care sector. Contact your Liko representative or visit [www.liko.com](http://www.liko.com) for further information and advice.



## Sit-to-stand

Sit-to-stand accessories can be mounted on Golvo without getting in the way during normal lifts. Easy installation; no tools required.



For more than 25 years, Liko has successfully specialized in the design, development and manufacture of safe, easy-to-use assistive devices: mobile lifts and overhead lifts, slings and accessories. Our specialization enables us to remain a step ahead. We gather knowledge and experience and apply it when we develop our products and methods.

Liko provides complete product service via our own service personnel and service partners. Please feel free to contact us for advice and guidance. Liko is a Swedish company. We are market leading in the Scandinavia and are represented by our own subsidiaries and distributors in more than 30 countries.

Liko is quality certified according to ISO 9001 and its equivalence for the medical device industry, ISO 13485. Liko is also certified according to environmental standard ISO 14001.



www.liko.com

**Manufactured by:**  
Liko AB  
Nedre vägen 100  
975 92 Luleå  
Sweden  
info@liko.se

**Liko Distributor in UK:**  
Liko (UK) Ltd • Brunel Centre  
Stroudwater Business Park, Stonehouse,  
Gloucester, GL10 3RU  
Tel: +44 (0)1453 823000  
Fax: +44 (0)1453 828844  
info@liko.co.uk

**National Distributor in USA and Canada:**  
Liko North America • 122 Grove Street  
Franklin, MA 02038  
Tel: +1 508-553-3993, 1-888-545 6671  
Fax: +1 508 -528-6642  
info@liko.com

